

TECHNICAL INFORMATION

OPTIMA EUROPA



REHABILITATION FOOT ORTHOSIS

Treatment of diabetic foot ulcers, post-operative/traumatic phases and foot fractures.

PRODUCT COMPONENTS

- 1- Leg
- 2- Lining
- 3- Single
- 4- Bottom
- 5- Velcro



The Optima line is indicated for the treatment of alterations in the morphology of the foot.

Optima Molliter S.r.l. declares that the Optima tutors are medical devices conceived and designed according to the state of the art, have been checked and made compliant with EU Regulation 2017/745 for the applicable parts. The braces comply with IWGDF's international guidelines for the management of diabetic foot ulcers. <https://iwgdfguidelines.org/guidelines/guidelines>.

THERAPEUTIC INDICATIONS

- Management of plantar, digital and interdigital ulcers
- Post-operative hallux valgus
- Post-operative hallux rigidus
- Post-operative hammer toes
- Post-operative care for Morton's neuroma
- Toe fractures
- Nail removal
- Plantar wart removal
- Post-operative forefoot
- Post-operative toes.
- Suitable for functional bandage support.
- Interdigital injury.
- Toe injury.

TECHNICAL FEATURES

Rehabilitation foot orthosis, thermo-mouldable heel counter, breathable and washable upper. Rigid anti-flexion and anti-torsion insole, Fully Rocking sole with movement control and pressure load redistribution.

TECHNICAL CHARACTERISTICS

MATERIALIALS:

UPPER: PU 70% + PES 20% + COT 10%

LINING POLYESTER: 100%

INSOLE INSERT: POLYAMIDE PA6 70% +GLASS FIBER 30%

BOTTOM SOLE: POLYURETHANE

PROTECTIVE SOCK: 90% COTTON, 7% POLYAMIDE, 3% ELASTANE

OPTIMA PRODUCTS ARE LATEX FREE.

CLOSURE: 2 VELCROS

PACKAGING: WEIGHT: GR 205 / MATERIAL: CARTON

OPERATIVE TEMPERATURE AND STORAGE: - 5° | + 45°

PRODUCTION LOT: AS DESCRIBED ON THE LABEL ON THE FRONT OF THE BOX

BASE COMPONENTS: MONOLAYER

ADDITIVES: INSOLE PLTM / OPTIMA KIT 3X3 PLANTAR INSOLE / PROTECTIVE SOCK

PRODUCT CLEANING


To clean the product, use warm or cold water and neutral soap; rinse once the operation is completed. Avoid any corrosive liquids such as alcohol, laundry detergents, bleach, ammonia, and solvents. Avoid using abrasive materials such as brushes, steel wool, and sand paper, or a knife blade to remove stains. Optima Molliter S.r.l. shall NOT be held liable for damages caused by abrasive products, by products which attack the materials, by direct exposure to heat sources.

Class I medical device according to Regulation (EU) 2017/745

In accordance with Recommendations 2-5-6-7c of the IWGDF International Guidelines on the Management of the Diabetic Foot. <https://iwgdfguidelines.org/guidelines/guidelines>

SIZES

	S 34-36	M 37-39	L 40-42	XL 43-45	XXL 46-48
A MM.	247	267	288	316	336
B MM.	101	108	115	119	131
C MM.	126	129	132	138	147
PESO GR. without insole	342	408	473	563	720

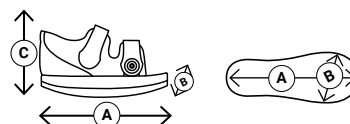

Optima Molliter Srl

Company Certified under: ISO 9001:2015 IT346378 - ISO 13485:2016 IT346377

Via Breda 19/21 62012 Civitanova Marche (MC) Italy T +39 0733 801060

info@optimamolliter.com / www.yda.it / www.optimamolliter.com

